



Knee Arthroscopy Surgery

Recommended accessories:
Donut head gel pad ×1
Torso section gel pad ×1
Arm rest ×2
Knee support



Rectal Surgery

Recommended accessories:
Arm rest ×2
Foot plate 1pair
Notched clamp ×2
Iliac support ×1



Orthopedic Procedures

Recommended accessories:
Donut head gel pad ×1
Torso section gel pad ×1
Orthopedics traction ×1



Chest/Abdominal Surgery

Recommended accessories:
Donut head gel pad ×1
Torso section gel pad ×1
Arm rest ×1
High position arm holder ×1
Universal dome shaped pad ×1



Orthopedic Procedures

Recommended accessories:
Orthopedics traction×1
Knee crutch for orthopaedic traction ×1
Installing rack of inserted rod ×1



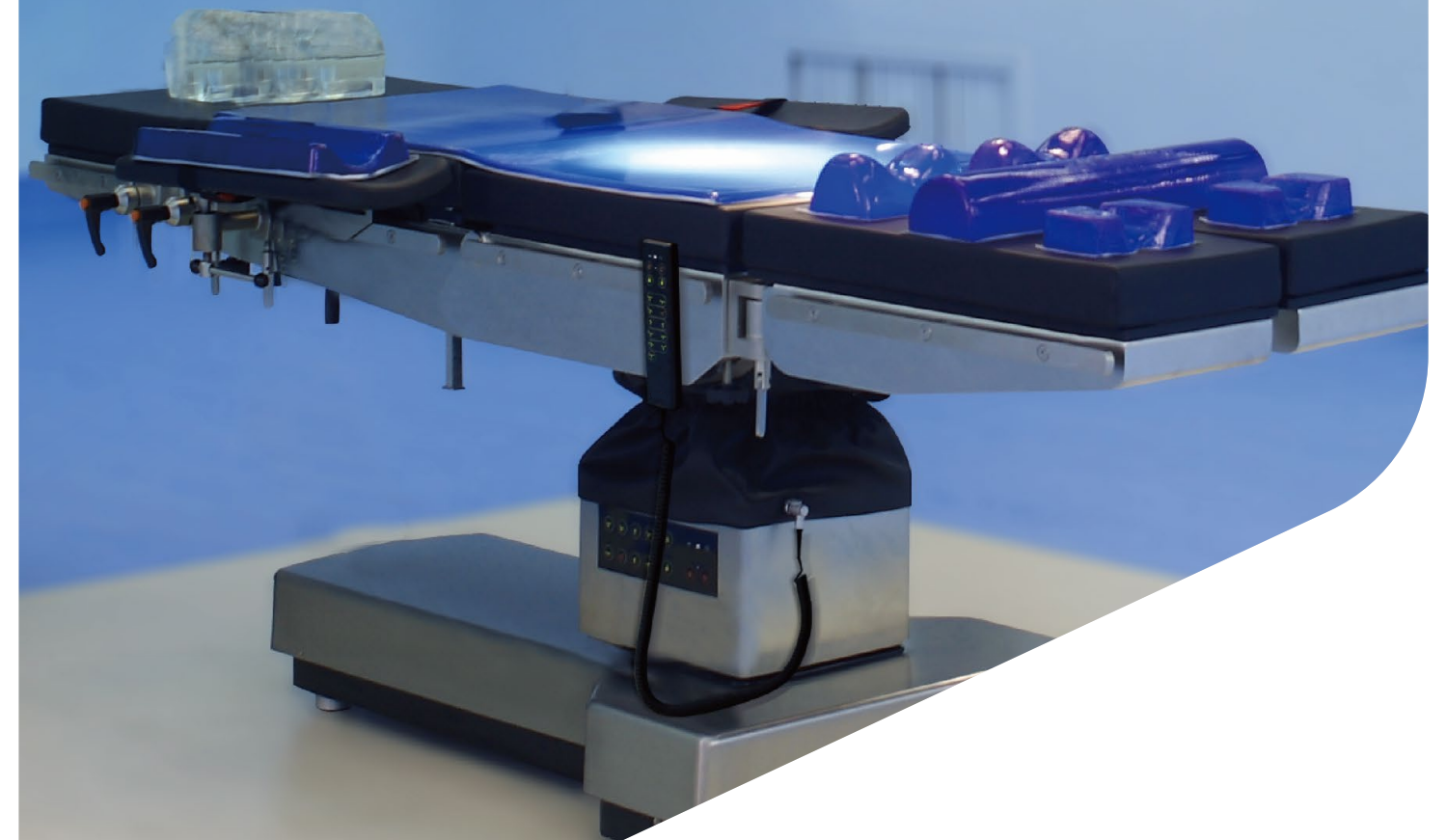
Orthopedic Procedures

Recommended accessories:
Orthopedics traction×1
Anti-tensile bracket for side-lying position×1
High position arm holder ×1

CE

Aegistab OP850

Electronic Operating Table



CORPORATE HEADQUARTER

HEYER Medical AG

Carl-Heyer-Str. 1/3, 56130 Bad Ems, Germany
www.heyermedical.de

HEYER Medical AG

Carl-Heyer-Str. 1/3, 56130 Bad Ems, Germany
Tel: +49 2603 791 0
Fax: +49 2603 791 200
info@heyarmed.com


HEYER Medical Co., Ltd.

No.10, Chaobai Street, Yanjiao Development Zone
065201, Sanhe City, Hebei, China
Tel: +86 10 5841 1198
Fax: +86 10 6371 8989
marketing@heyarmed.com

AOP850-2022.05



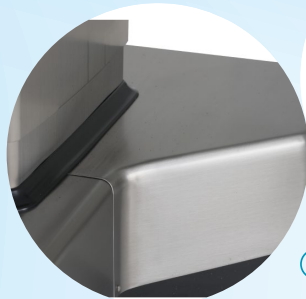A high and multi-performance electronic operating table with a maximum load of 350kg, meeting the requirements of different surgeries.



① Premium Quality Material

- All exposed parts are made from stainless steel SUS304
- Easy to clean and disinfect



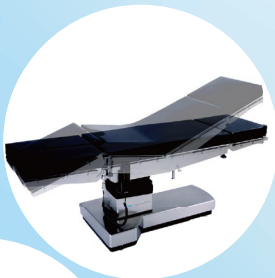
② Built-in Elevator Kidney-bridge

- Has no negative effect on X-ray radiolucency and is able to support kidney and gall bladder surgeries efficiently




③ Perfect for X-Ray Imaging

- Head and leg sections with excellent radiolucent table guarantee high quality X-ray imaging
- Exchangeable head & leg sections
- Optional with 300 mm longitudinal shifting



④ Surgical Configuration

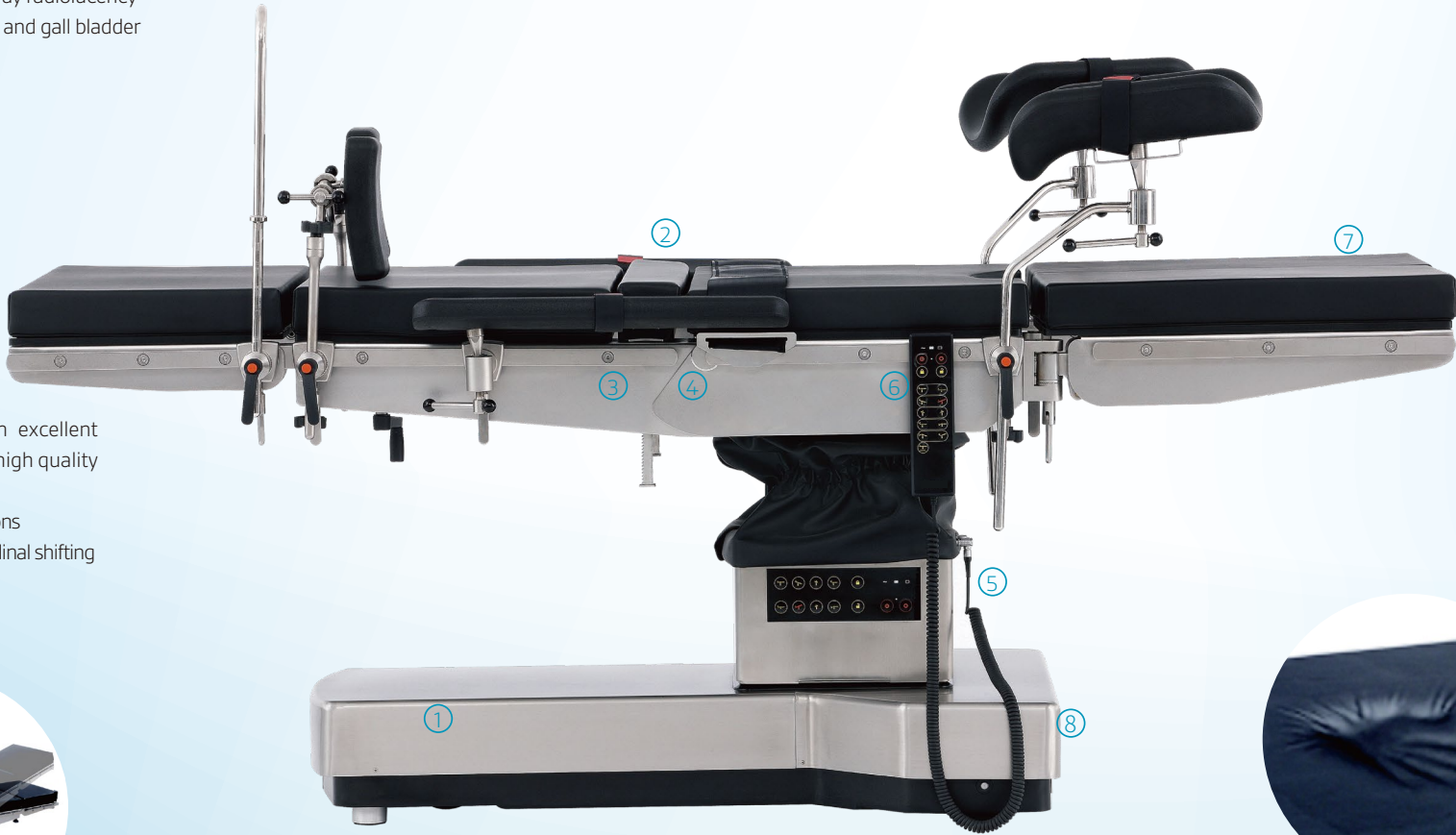- Has one-key reset, flex and reflex function
- Table angles and section lengths can be adjusted within a wide range with 6 standard and various optional accessories, meeting the need of various surgeries



⑤ Ergonomic Adjustment

- Convenient design for positioning of leg plate and head plate, able to control leg and head plates with one hand
- Imported gas springs guarantee smooth movement
- Imported round-hole clamps enable easy lock and unlock accessories, ensuring safety during operation
- Adjustable table height (650-1,000mm) and T-shape base provide comfortable working environment



⑥ Safety

- Auto-lock system: table base will be locked automatically once the position was adjusted by controller, preventing possible misoperation
- Two control modes: manual and column control panel, allowing for backup if one is broken, ensuring the safety of surgery
- Control panel will power off automatically when there is no operation over 1 minute, avoiding any possible misoperation



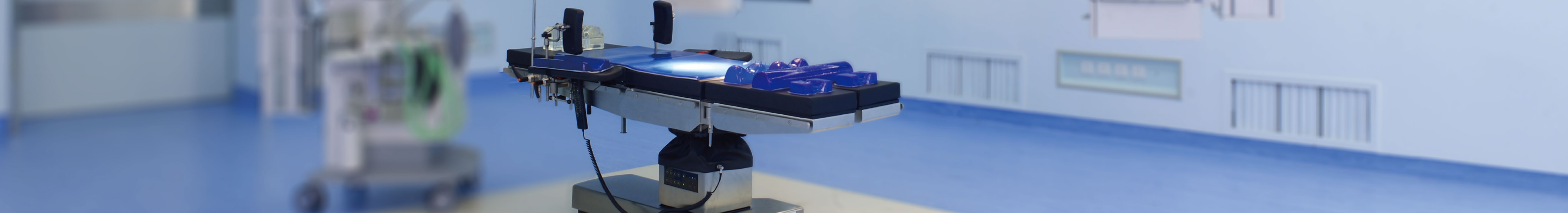
⑦ Pads

- Have soft texture and stuffed with memory foam
- Effectively reduce possibility of patient bedsores




⑧ High Efficiency Built-in Battery

- Imported from Panasonic
- Maximum 15 days working time
- Maximum 85 normal operations duration
- No request for AC power to increase the safety of surgery



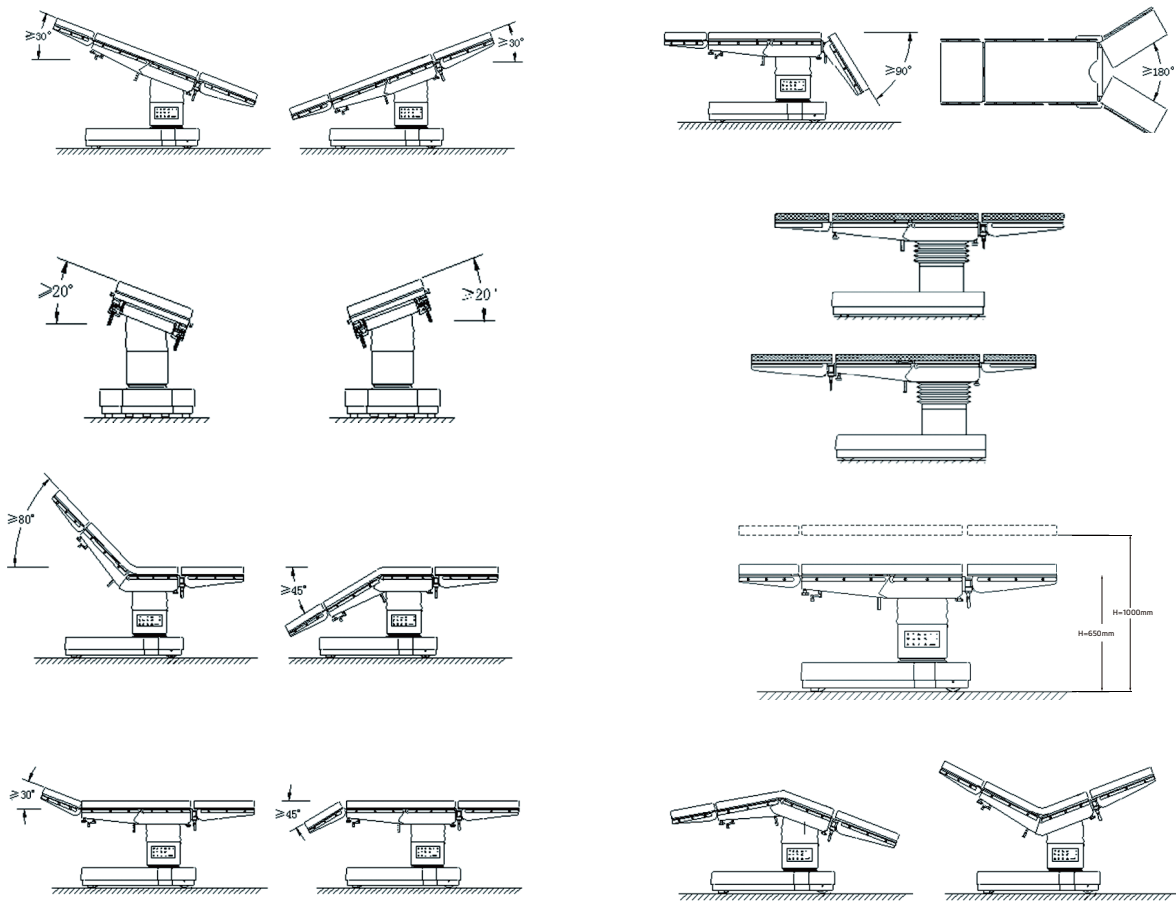Wide adjustment range for all positions, satisfying the requirements for different surgeries.



OP850 is a high quality hydroelectric operating table which satisfies european clinical safety requirements for convenient and flexible operation system. It will be your best choice for modern operation room.


Technical Specifications

| | |
|-------------------------------------|------------------|
| Table Length(mm) | 2,080 |
| Table Width(without side rails)(mm) | 520 |
| Table Height(without pads)(mm) | 650 to 1,000 |
| Trendelenburg / Reverse | 30°/30° |
| Lateral Tilt (left/right) | 20°/20° |
| Back Plate Up/Down | 80°/45° |
| Head Plate Up/Down | 30°/45° |
| Leg Plates Down | 90° |
| Leg Plates Split(°) | 180 |
| Flex/Reflex(°) | 225 / 100 |
| Kidney-bridge Elevator Up(mm) | 128 |
| Power Input(Hz) | AC 90V-264V50/60 |
| Cell Output(AH) | DC 24V12 |
| Weight Capacity(kg) | 350 |

Optional Accessories

Features:

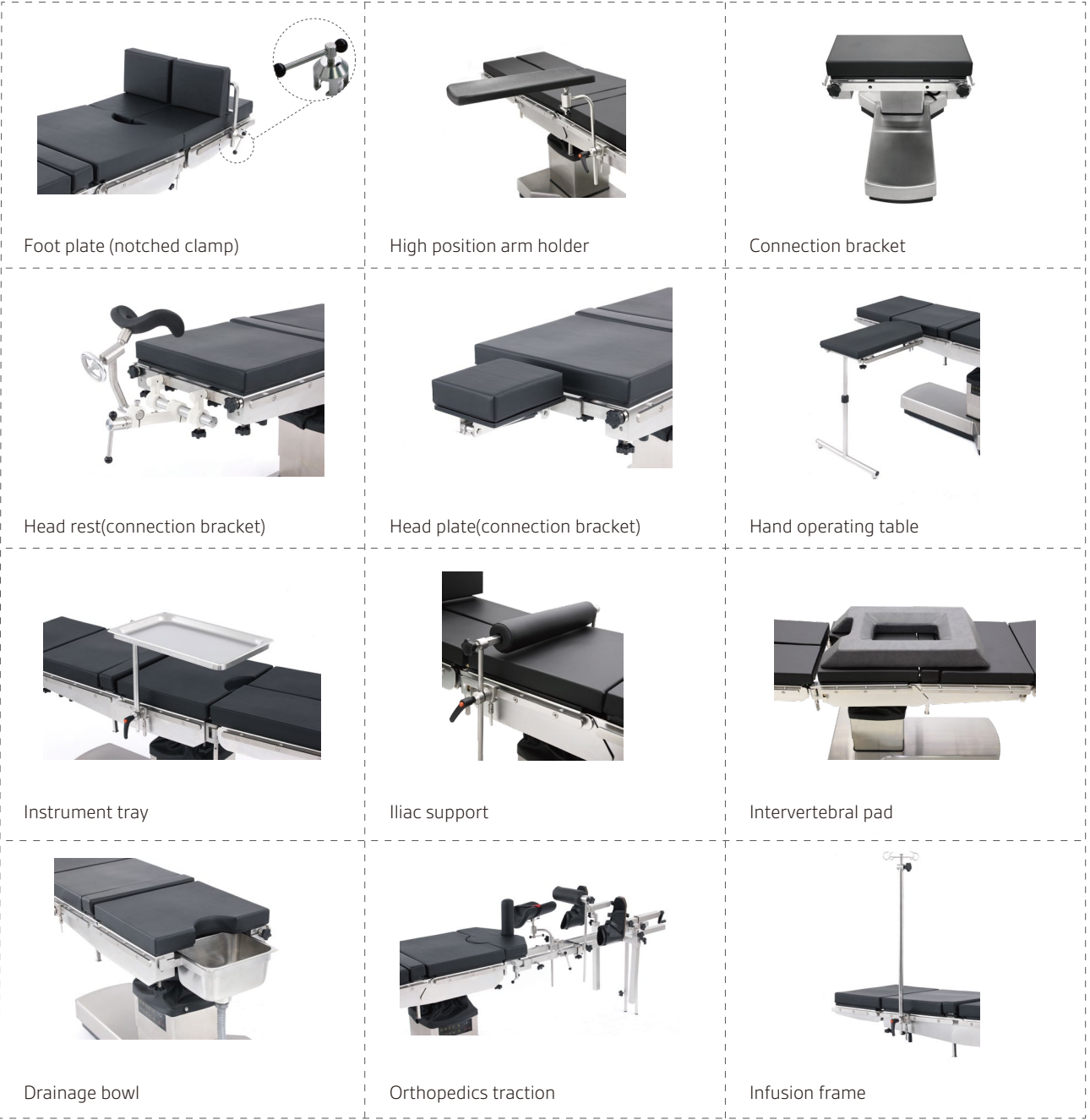All table accessories are made of SUS304 stainless steel. All the surfaces are easy to clean and disinfect. With 6 standard and 13 optional accessories, the surgical bed can support a multitude of body positions and procedures, meeting the requirements of most surgeries.



Orthopedics Features





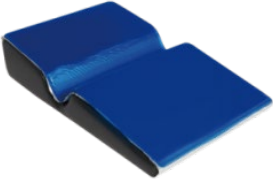



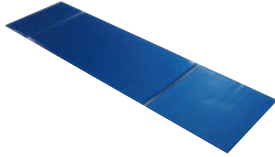
The orthopedics traction is designed based on EU's requirements, meeting the requirements of most orthopedics surgeries. All parts are made of SUS304 stainless steel with high strength. Therefore, the surgical table is very stable during surgeries and it is also X-ray and C-arm radiolucent.

Standard Accessories

| | | |
|---|---|---|
|  |  |  |
| Round-hole clamp | Anesthesia screen | Body strap |
|  |  |  |
| Arm rest | Lateral support | Knee crutch |

Gel Pad
Functions & Features:

The gel pad which made from viscoelastic polymer gel, offers stability and comfort during positioning. It has uniform dispersion of pressure which can avoid bedsores. It is easy to clean and X-ray radiolucent. It can support surgery with working environment temperature from -10°C to 50°C. It is also biocompatible with human skin and tissues. All these features make the gel pads cost-effective and durable.

| | | |
|---|---|---|
|  |  |  |
| Prostrate head pad | Donut head pad | Horseshoe head pad |
|  |  |  |
| Ophthalmic head pad | Lateral pad | Universal dome shaped pads |
|  |  |  |
| Contoured arm/Leg pad | Heel pad | Torso section |

Surgery Positions:



Chest/Abdominal Surgery
Recommended accessories:
Donut head gel pad ×1
Torso section gel pad ×1
Foot plate 1 pair × 1
Notched clamp × 2
Arm rest × 2



E.N.T Surgery
Recommended accessories:
Donut head gel pad ×1



Urology/OB/Gynecology
Procedures Surgery
Recommended accessories:
Donut head gel pad ×1
Torso section gel pad ×1
Knee crutch × 2
Drainage bowl × 1
Arm rest × 2



Spinal Surgery
Recommended accessories:
Prostrate head gel pad × 1
Intervertebral pad × 1
Arm rest × 2
Torso section gel pad ×1
Universal dome shaped pad ×1



Shoulder Operation
Recommended accessories:
Arm rest × 2
Head rest × 1
Connection bracket × 1



Neurosurgery
Recommended accessories:
Arm rest × 2
Head rest × 1
Connection bracket × 1
Torso section gel pad ×1
Universal dome shaped pad ×1